



producing quality



microscope

# OP-C16

# OP-C16



## ERGONOMICS AND VERSATILITY

The OP-C16 microscope manufactured by OPTOMIC fulfils the most demanding requirements.

The OP-C16 microscope offers one of the widest ranges of lighting systems and beam splitters, as well as image and video capture and treatment devices in the market.


Its excellent optical and mechanical features, together with its great mobility and versatility, make the OP-C16 microscope the ideal equipment for any treatment.

The sturdiness and stability of the OP-C16 microscope, as well as the freedom of movement, enable the head to be placed in any desired position to allow the specialist to adopt an ergonomic posture.

## COMPACT AND VERSATILE

Its arm can be folded to minimum dimensions, thus resulting very practical in the specialist's office as it occupies very little space.

## VERSATILITY FOR ALL TECHNIQUES

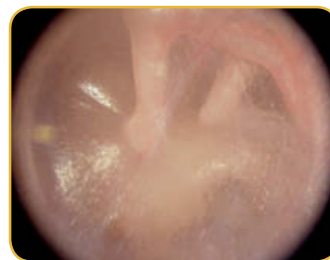

OPTOMIC thanks Dr. Luis García-Ibañez for contributing these photographs.

The OP-C16 microscope is designed and built using the most modern and precise manufacturing systems.

Soft Move is our exclusive movement and braking system, which provides OPTOMIC equipment with a play-free turning mechanism and an excellent mobility, smooth and precise. Soft Move is supported by the great sturdiness and stability of the highly resistant conical cylinder bearings and a unique braking

system. The large central knob of this braking system acts by freeing or braking the turning movements with absolute precision, enabling even a total blockage.

The Soft Move system provides the specialist with extremely smooth and accurate movements, permitting the doctor to adopt a comfortable and ergonomic position during examination and surgery.

The head, with its swinging and directional movement, has a brake in each position, as well as a micrometric focusing knob. The magnifications changer is mounted on angular contact bearings, providing free lubrication maintenance, inexistence of play and an absolute accuracy in the rotation and vision of the equipment. This system is exclusive of OPTOMIC equipment.

## supports

### SUPPORTS

#### Rolling Floor Stand

The floor stand is of such reduced dimensions that the microscope can be positioned close to the ENT equipment without interfering with the rest of the apparatus. This very stable rolling floor stand has 5 double rubber wheels with independent brakes. These features give our equipment a great mobility as well as perfect balance and safety.

#### Fixed Floor Mount

This floor mount confers the microscope with a great stability. Its base has a diameter of only 360mm and is directly fixed to the floor, making it possible to position it very close to the operative unit.

#### Wall Mount

With the long arm of this space-saving design, the specialist is able to work at a considerable distance from the wall. When not in use, it can be folded against the wall. Using an optional extension, a maximum 2.10 m radius can be achieved, thus guaranteeing the desired reach.

#### Ceiling Mount

This type of mount provides comfort and agility in the specialist's office and surgery rooms.



# INTELLIGENT COMPATIBILITY AND INNOVATION

## Beam Splitters

OPTOMIC offers a complete range of beam splitters to meet the needs of every specialist.

### OP-TV6

It includes a compact CCD camera with an independent focusing system.



### OP-TV8 / OP-TV10 HD

Compact system for video-microscopy, with integrated high resolution or HD camera. It incorporates white balance control and power switch on the beam splitter.



### OP-TV DC + DIGITAL CAMERA

With an exclusive system to adapt it to digital photographic cameras.



### OP-TV CLICK

It adapts to any endoscopy CCD camera to obtain microscopic images with digital output of outstanding resolution.



## lighting LIGHTING

The OP-C12 microscope has three different illumination options, to adapt to any needs. All light sources have a connector for fiber optic cable for endoscopes, becoming a multidisciplinary device.

### FIBROLUX 150

Double easy-to-change 150w halogen lamp. Lamp protection system to lengthen its life.



### FIBROLUX LED HP

The latest technology in LED illumination. High lighting efficiency. It has an additional filter for a warmer light, as well as a remote switch and power outputs for accessories such as monitor, camera or laptop.



### FIBROLUX X100

Bright illumination and outstanding colour reproduction provided by the 100w Xenon lamp. Lamp protection system.





## Binoculars

The OP-C16 microscope offers three options:

- Wide field straight binocular (standard)
- 45° inclined binocular (optional)
- 0-240° inclinable binocular (optional)



## Objectives

- Standard 200mm.
- Optional 250mm, 300mm and 400mm.

Thanks to its large diameter lenses, it favors the coaxial illumination of the observation axis, offering a homogeneous and bright illumination field, as well as a better light transmission and an optimum lighting of the area explored.



## Eyepieces

Standard: 12.5x widefield eyepieces.

Optional: 10x and 16x widefield eyepieces.

Large diameter lenses.

The optical coating offers a high resolution, contrast and definition, providing bright and chromatically real images, and an extraordinary viewing quality.



## Effective Magnifications

5-step Galilei changer, permitting magnification from 1.0x up to 34x, depending on the chosen combination of objective, binocular and eyepieces.

OPTOMIC designs and manufactures all of its optical components with the highest quality and precision standards, using the best glass types and coatings, resulting in equipment of great reliability and excellent optical quality.

# digital microscopy



## DIGITAL MICROSCOPY

### Digitally Soft

Digital microscopy integrates the features of the OP-C16 microscope and the advantages of the digital image, video and data capturing system, Digitally Soft.



Our Digitally Soft image processing software makes it possible to obtain reports, export images and videos and compare the pre- and post-treatment conditions. It also admits statistical studies, area calculations, compares measurements and has other functions which permit the patient's data to be processed in a simple and agile manner.

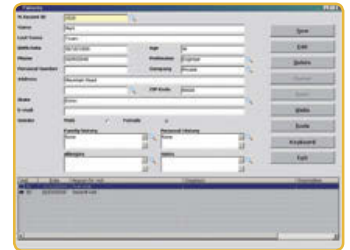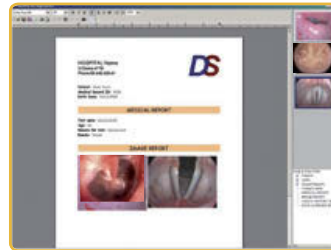

# Digitally Soft

## Digitally Soft



# microscopy-endoscopy

## MICROSCOPY-ENDOSCOPY


### Two in one

In the OPTOMIC digital video-microscope we find the integration of the techniques of microscopy and endoscopy thanks to the compatibility of the lighting and image systems.

Simply connecting a fiber optic cable to the output of the microscope's light source and using the microscope's camera (CCD 1428/714) will result in perfect endoscopic equipment, optimizing both space and resources.



## DIMENSIONS

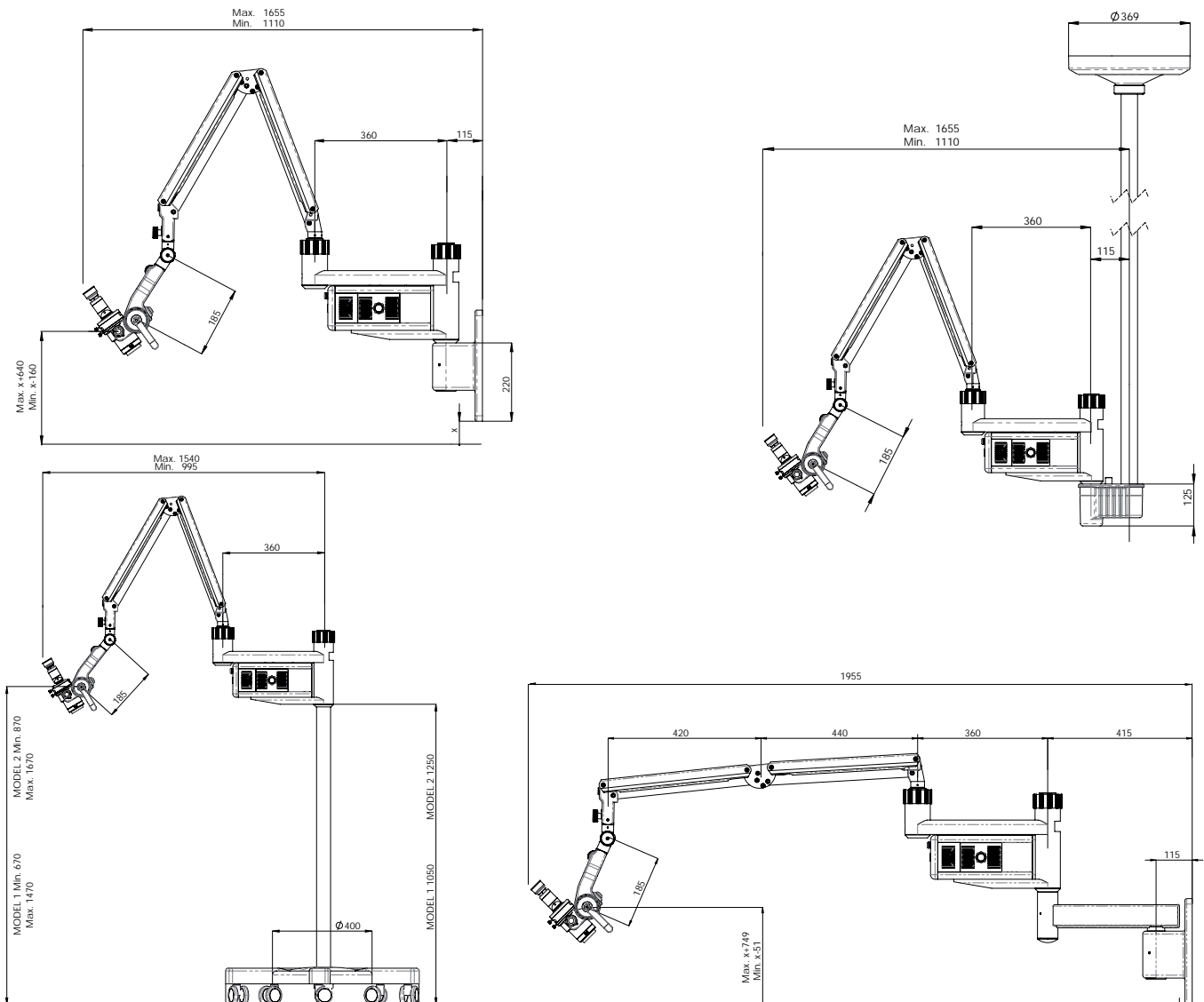


# technical features

## TECHNICAL FEATURES

OP-C16

HEAD	Stereoscopic head with micrometric focus
BINOCULAR	Wide field straight binocular (standard) 45° inclined binocular (optional) 0-240° inclinable binocular (optional)
ILLUMINATION SYSTEM	Halogen, LED or Xenon light source
LIGHT CONDUCTION	Fiber optic cable
EYEPIECES	12.5x widefield eyepieces. 10x and 16x optional
MAGNIFICATIONS CHANGER	5-steps Galilei
OBJECTIVE	200mm standard objective. 250mm, 300mm and 400mm (optional)
HEAD AND ARM	With brakes and movements using angular contact bearings and high resistance conical cylinder bearings, Soft Move system
VOLTAGE	100/240V 50/60Hz
MOUNTING OPTIONS	Floor, wall or ceiling mount
ACCESSORIES	Beam splitters and video systems Digitally Soft software Monitor support Laptop tray Camera tray

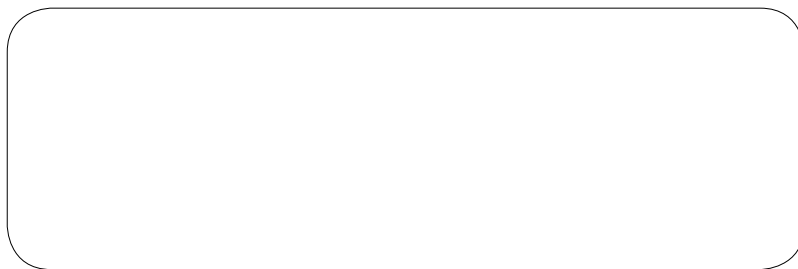

Specifications are subject to change without notice

